

All-in-one design  
comes in pairs



## **Compact** Hysteroscopes

Slim, compact design  
in 5.0 mm and 3.8 mm diameters

# Compact Hysteroscopes

Slim, compact design  
in 5.0 mm and 3.8 mm diameters

The slim, all-in-one design of our compact hysteroscopes from Richard Wolf does not only offer aesthetic appeal, but also practical benefits.

Good handling, optimized reprocessing and a long usage life meet internal qualities such as a brilliant image, true continuous flow and a large, straight working channel.

Challenges are easily overcome with two sizes of compact hysteroscopes:


## All-in-one comes in pairs

### Well-protected


The durable, compact construction protects surely and steadily the premium-quality PANOVIEW PLUS optics. This increases the lifespan greatly compared to dismantling instruments and ensures a sustained, optimal image quality.

### Lots of room with minimal impact


Even with a relatively thin outer diameter, both models are equipped with a large 5 Fr. working channel, which allows the passage of various instrumentation.

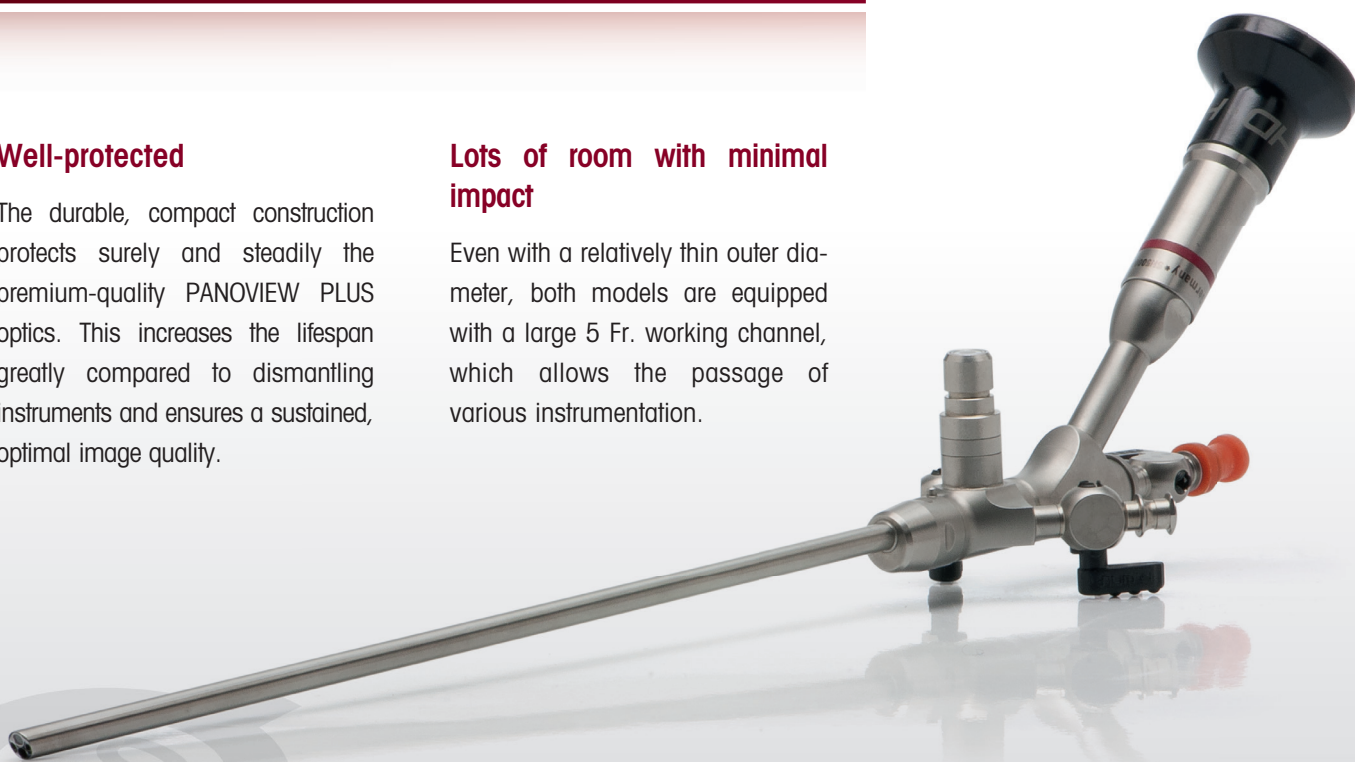

3.8 mm

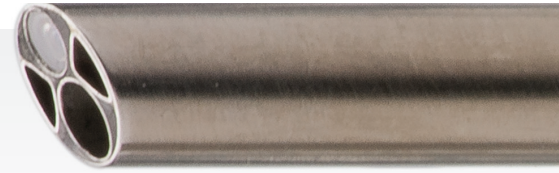

5.0 mm

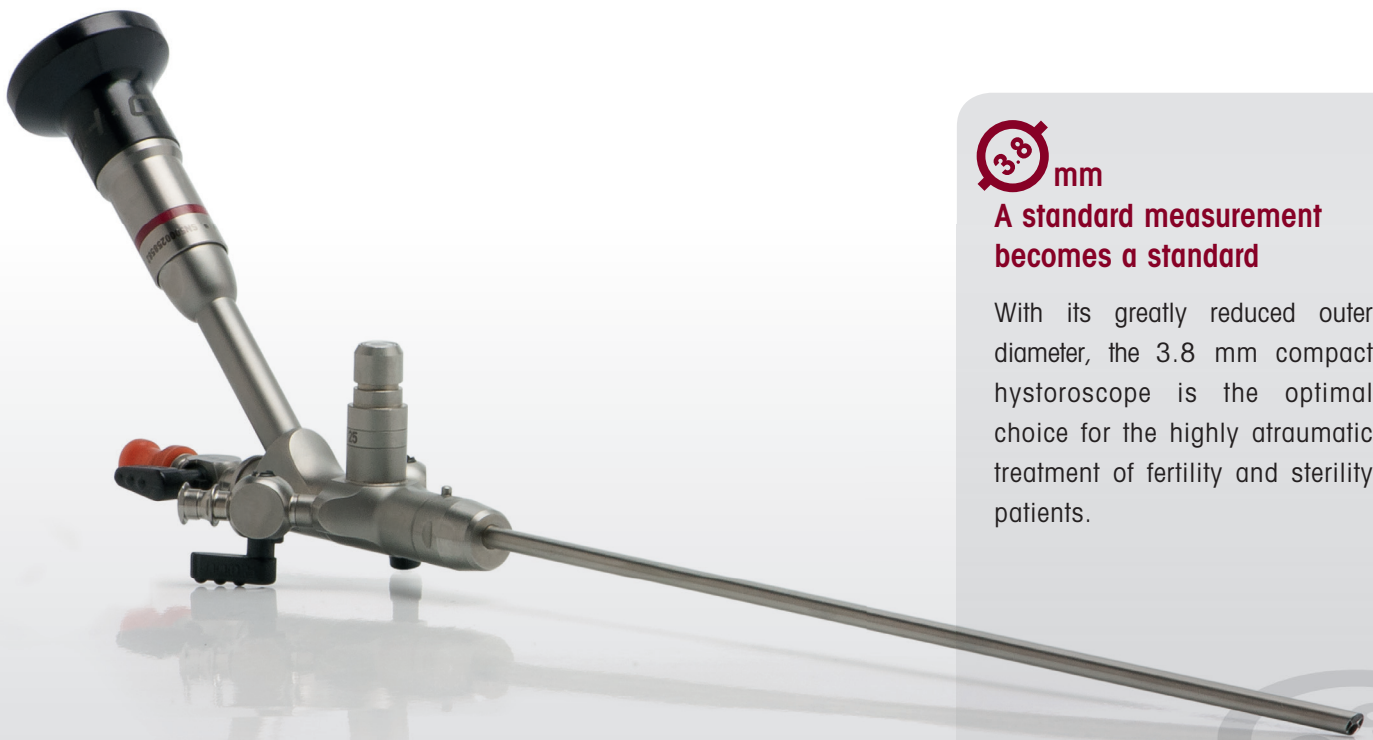


**A standard measurement becomes a standard**

With its greatly reduced outer diameter, the 3.8 mm compact hysteroscope is the optimal choice for the highly atraumatic treatment of fertility and sterility patients.

**True continuous flow**

Two separate in and outflow channels ensure a constantly clear sight of the operating field - also when instruments are passed through the working channel.

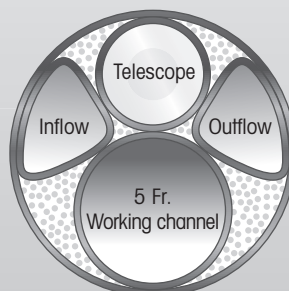

# Compact Hysteroscopes

Slim, compact design  
in 5.0 mm and 3.8 mm diameters

## Simply straight on

The straight 5 Fr. working channel ensures hand instruments remain protected from bending during operation. The attachable ergonomic handle allows comfortable and easy handling of the endoscope.



Compact hysteroscope	8986.401	89864021
Outer diameter	5.0 mm round	3.8 mm round
Working channel	5.0 Fr.	5.0 Fr.
Direction of view	20°	30°
Length	348 mm	352 mm
Working length	217 mm	217 mm
Weight	158 g	123 g

### Compact hysteroscope 5.0 mm

with continuous flow, direction of view 20°, straight working channel, for instruments up to 5 Fr. .... 8986.401

### Compact hysteroscope 3.8 mm

with continuous flow, direction of view 30°, straight working channel, for instruments up to 5 Fr. .... 89864021

### Accessories

<b>Ergonomic handle</b> .....	8986.251
<b>Biopsy forceps (Spoon)</b> .....	8642.6002
<b>Grasping forceps with alligator jaws</b> .....	8642.6502
<b>Micro scissors</b> .....	8642.6812
<b>Biopsy punch</b> .....	8642.6312
<b>"Bipotrode" electrode, bipolar 5 Fr.</b> .....	823005
<b>Reprocessing basket</b> .....	38044111
<b>Fiber light cable set, 2.3 m long, fiber bundle ø 2.5 mm</b> .....	806625231
<b>Rubber cap "RIWO", Red-orange, pack of 10</b> .....	88.01
<b>Stopcock cpl., Pack of 5</b> .....	896.0002