






CANTEL IS PROUD TO OFFER THE SURESTORE™ ENDOSCOPE STORAGE & TRANSPORTATION

British designed and manufactured by our experienced engineers and technical staff, the SURESTORE[™] Endoscope Storage and Transportation solution offers cutting edge long-term 'active' aseptic storage and transportation of endoscopes (including GI and theatre-based scopes). Conditioned directly from an endoscope washer disinfector such as RAPIDAER[™] Endoscope Reprocessor, and completely negating the need for time consuming air drying, the SURESTORE[™] Endoscope Storage and Transportation solution validates each scope for storage of up to 100 days.

The unique design and technology of the SURESTORE[™] Endoscope Storage and Transportation solution provides a 'patient-ready' packaged endoscope within 10 minutes.

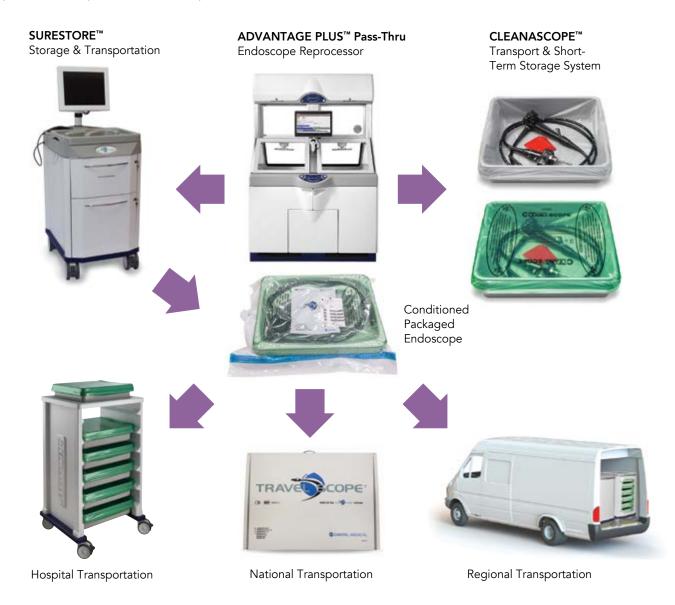



THE COMPLETE CIRCLE OF PROTECTION

As the global vanguard in infection prevention, **only Cantel delivers the Complete Circle of Protection**, a full-value, proactive partnership dedicated to helping you remove risk, streamline operational efficiencies and optimise your success.

DRY & STORE Bacterial growth poses significant risk to endoscopes during storage. Our Transport, Drying and Storage solutions are designed to protect valuable inventory, reduce cross-contamination touch points, eliminate moisture in the endoscope channels and control humidity. Higher humidity or moisture in endoscope channels is known to aid bacterial growth.

ENDOSCOPE STORAGE & TRANSPORTATION



SURESTORE[™] Endoscope Storage and Transportation solution enables scopes to be easily prepared in advance for emergency situations, out of hours, weekend, and off-site endoscopy work. The scopes can be transported to various departments within a hospital, or to secondary sites, using the TRAVELSCOPE[™] Transportation Case, or the TRANSCOPE[™] Lockable Carts and STORASCOPE[™] Lockable Carts.

- Validated storage for 100 days
- Independent monitoring system
- Reduced scope wear and tear
- Flexible transportation solution
- Fully traceable cycle data
- Quick and simple to use
- Secure and robust storage
- Space saving storage solution
- Independently validated by Biotech Germande
- Tested with multiple endoscope brands



Model

SURESTORE[™] Endoscope Storage & Transportation

Conditioning time (start to finish)

Approximately 8 minutes (lumen scopes) Approximately 2 minutes (non-lumen scopes)

Dimensions

545mm (w) x 690mm (d) x 1220mm (h) (960mm exc screen)

Power requirements 13 Amp 220 - 240V

Operating requirements

Completely mobile on 4 independently lockable castors

Installation space

Unit does not require any additional ventilation space and can therefore be positioned against a wall or between other equipment

Validation

Biotech Germande and NFS 98-030

ORDERING INFORMATION		
PART NUMBER	DESCRIPTION	UNITS PER BOX
100600UK/DE/F/NL/IT (landcodes depend on destination kit country)	1 X SURESTORE [™] Storage/Transportation Unit 12 X SURESTORE [™] Supportive trays 1 X 5 Litre N-SURE [™] Solution 50 X SURESTORE [™] Outer Pouches 45 X SURESTORE [™] Sterile Collection Pads 50 X CLEANASCOPE [™] Liners	1
100602	N-SURE [™] Solution 5L	1
100604	SURESTORE [™] Outer Pouches	50
100728	SURESTORE [™] Sterile Collection Pads	225
100606	SURESTORE [™] Supportive Tray	1
700344	SURESTORE [™] Tracking Labels	400
200018	SURESTORE [™] HEPA Air Filter	1
103580	$CLEANASCOPE^{\texttt{M}} \text{ Sterile Liner Set (red and green covers with base liner)}$	200
100521	TRAVELSCOPE [™] Transit Case	1
103551	TRANSCOPE [™] Lockable Cart 6 Tray (contoured trays)	1
103543	TRANSCOPE [™] Lockable Cart 6 Tray (standard trays)	1
103544	TRANSCOPE [™] Lockable Cart 6 Tray (standard trays with lids)	1
103553	STORASCOPE™ Lockable Cart 10 Tray (contoured trays)	1
103547	STORASCOPE [™] Lockable Cart 10 Tray (standard trays)	1
103548	STORASCOPE™ Lockable Cart 10 Tray (standard trays with lids)	1

 $\label{eq:random} RAPIDAER^{\texttt{W}}, SURESTORE^{\texttt{W}}, CLEANASCOPE^{\texttt{W}}, TRAVELSCOPE^{\texttt{W}}, TRANSCOPE^{\texttt{W}}, STORASCOPE^{\texttt{W}} and N-SURE^{\texttt{W}} are trademarks of Cantel (UK) Limited$

© 2017-2021 Cantel. All Rights Reserved. The Complete Circle of Protection diagram is a trademark owned by Cantel Medical Corp. and is subject to copyright protection owned by Cantel Medical Corp.



Medivators BV | t: +31 (0)45 547 1471 | e: csdistributors@medivators.com

