

## PANETTI ENDOSCOPIC / MICROSCOPIC INSTRUMENT SET - VISUAL AID

Suction dissector for ear drum, left side, 3 mm    Suction dissector for ear drum, right side, 3 mm    Suction dissector for ear drum, left side, extra fine    Suction dissector for ear drum, right side, extra fine    Suction dissector for sinus tympani, left side, 6 mm    Suction dissector for sinus tympani, right side, 6 mm    Suction dissector for atticus, left side, 6 mm    Suction dissector for atticus, right side, 6 mm    Suction micro hook, 1.0 mm    Suction micro hook, 1.5 mm    Suction micro hook, 2.0 mm



Dissection and suction of the cholesteatoma around the ear drum

Dissection and suction of the cholesteatoma in the sinus tympani

Dissection and suction of the cholesteatoma in the atticus

Cutting through the incudo-stapedial-joint; lifting of the footplate

Suction separator suction inner side 2.5 mm    Suction separator suction outer side 2.0 mm    Suction separator suction outer side 2.5 mm    Suction knife suction left side    Suction knife suction right side    Suction sickle knife    Suction curette    Suction dissector left curved, blunt    Suction dissector right curved, blunt    Preyer suction probe 90° angled, blunt



Dissolving and removing of the ear drum

Incision of the external meatus; preparation of the meatal flap

Cutting tissue

Removal of bone structure

Removal of cholesteatoma in the antrum

Exploring the antrum

**New instruments – New possibilities**



10-815-12 Suction sickle knife for ear drum, 45° angled, 12,5cm



10-816-12 Suction curette, angled, 12,5 cm



**New instruments – New possibilities**

10-817-12



Suction dissector,  
upcurved left,  
blunt

10-818-12



Suction dissector,  
upcurved right,  
blunt

10-819-12




Preyer suction finder,  
blunt



### Panetti ear drum suction dissector

- Very delicate shaft
- Interesting for small areas or children
- Better view around the corner

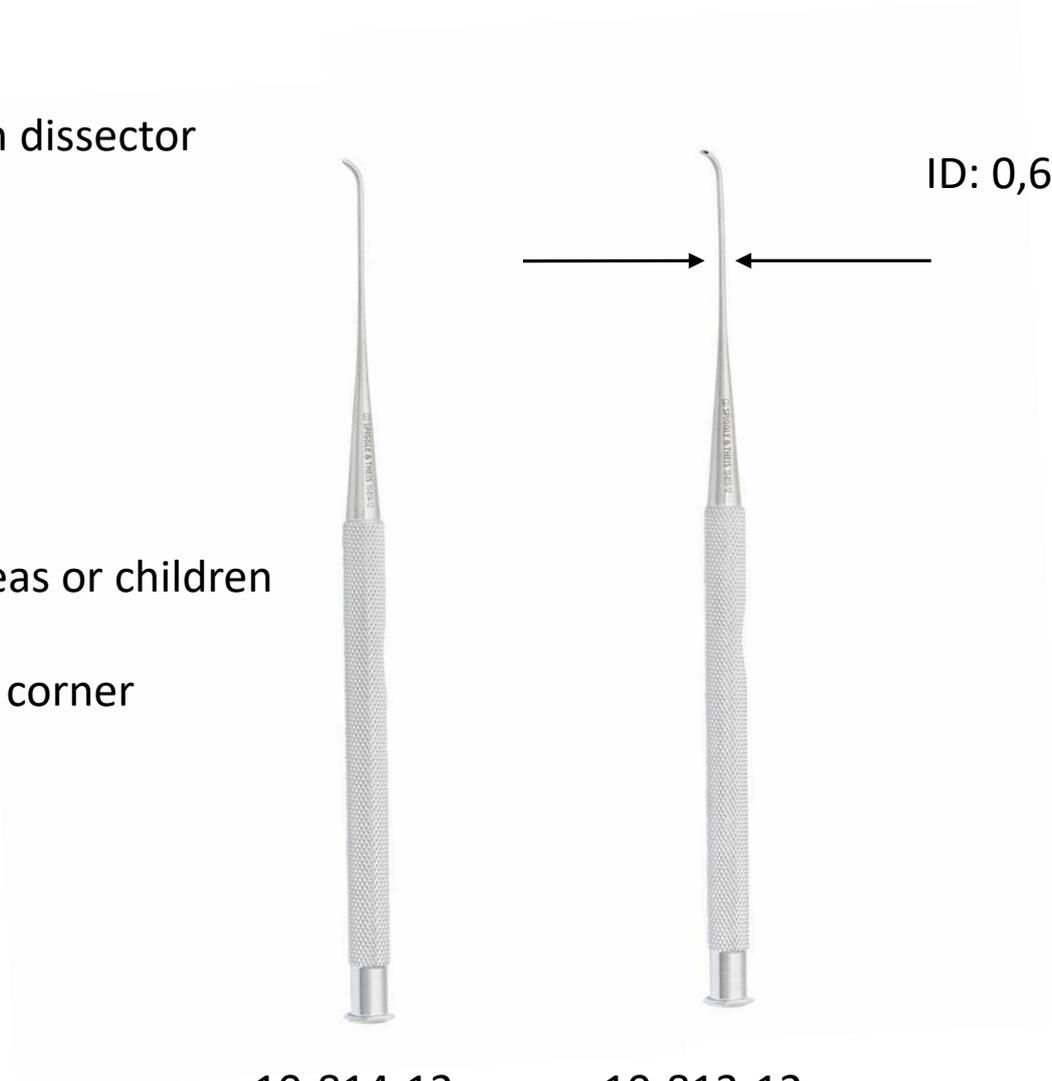

10-814-12

10-813-12